LANAFORM®

Dolce
Beauty


Exfoliating action

Action exfoliante

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
SCHEME A


1

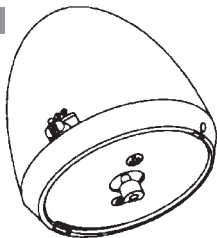
2

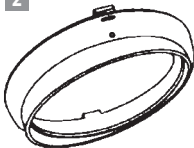3

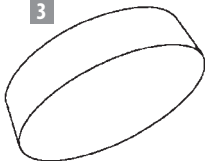
4

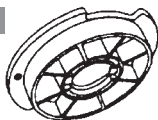5






6

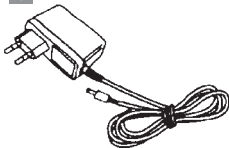7

SCHEME B

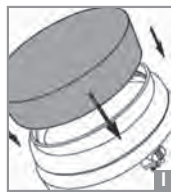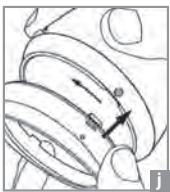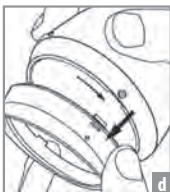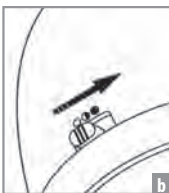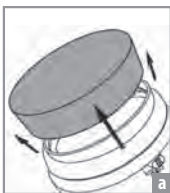EN - INTRODUCTION

Thank you for choosing Dolce by Lanaform®.

Providing real foot care, this Dolce pedicure device will make your feet naturally soft. Ideal for damaged, dry and rough feet, this exfoliating care has been designed to smooth the skin of your feet and heels. Thanks to its rotating abrasive disk, dead skin on your feet will be easily eliminated.

PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR DEVICE, ESPECIALLY THIS BASIC SAFETY INFORMATION:

- This device is not intended for use by people, including children, with limited physical, sensory or mental abilities, or by people lacking experience or knowledge, unless they have been supervised or instructed beforehand on how to use the device by someone responsible for their safety. Children should be supervised to ensure that they do not play with the device.
- Only use this device according to the instructions given in this manual.
- Check that your mains voltage is the same as that of your device.
- Keep the device away from heat sources such as stoves, radiators, etc.
- Do not leave the device in direct sunlight.
- Do not connect or disconnect the plug or electrical transformer of the device with wet hands.
- If the power cable is damaged, it must be replaced by a similar cord available from the supplier or the latter's after-sales department.
- Unplug the device before cleaning it.
- Do not wash the entire unit with water and do not immerse it in water. Consult the cleaning instructions in this manual.
- Unplug a device that has fallen into water immediately before retrieving it.
- An electrical device should never remain connected unsupervised. Disconnect it when you are not using it.
- Keep the power cable away from hot surfaces.
- Never use this device in a room where aerosol products (sprays) are used or in a room where oxygen is being administered.
- Do not use this device under a cover or under a cushion. Excessive heat can cause fire, electrocution or injury.
- Never use this device in a room whose temperature exceeds 40°C.
- If the device is placed in a room where the temperature is low, you are advised to wait for several minutes before using the device in a warmer environment. Indeed, drops of water may appear due to condensation, in which case your device may not work properly.
- Do not put items of clothing or towels under or on top of the device.
- This device is designed for domestic use only.
- If you are concerned about your health or are taking medical treatment, consult a doctor before using this device.
- If you experience any pain when using this device, stop using it immediately and consult your doctor.
- Do not use this device for more than 30 minutes at a time. Wait for one hour between each use.
- Never use this device on a body part that is swollen or inflamed or if you have a skin rash.
- Do not use this equipment without consulting your doctor in the following cases: if you are pregnant or breast-feeding, if you have diabetes, if your

skin has tattoos or moles, if you have skin cancer, if you have unusual skin discolouration, if you take medicine or if you have a foot problem such as nail fungus or over-sensitive skin. If you are unsure, consult your doctor before using Dolce.

COMPONENTS OF THE DEVICE (SEE ILLUSTRATION «A» AT THE BEGINNING OF THESE INSTRUCTIONS):

1. Main body (1)
2. Protective ring (2)
3. Cover (3)
4. Receptacle for callused skin (4)
5. Abrasive disk (5)
6. Mains adaptor (6)
7. Cleaning brush (7)

INSTRUCTIONS FOR USE (SEE ILLUSTRATION «B» AT THE BEGINNING OF THESE INSTRUCTIONS)

HOW TO USE

For an optimum result, soak your feet in warm water for 15 minutes and dry them before using Dolce. Rinse them to finish the treatment and moisturise them with a nourishing cream.

1. Open the cover. (a)
2. Connect the adaptor to the electric socket and connect the other side to the socket on the main body.
3. Put the lever in position 1. (b)
4. Take the main body and apply it onto your feet using circular movements across your calluses to remove dead skin. Do not stay in the same place but rather move the device across a large surface area of skin so as to not irritate it.
5. If you have thick calluses, you can put the lever in position 2. (b)
6. Keep the upper part of the protective ring aimed at your face and adjust the angle of the main body's nozzle to keep residues from escaping towards your face.
7. After use, put the lever in off position and unplug the adaptor from the electric socket and device.

DISMANTLING AND ASSEMBLY

1. Turn the protective ring clockwise to separate it from the main body. (c & d)
2. Lift the abrasive disk vertically to separate it from the receptacle for callused skin. (e)
3. Vertically remove the receptacle to separate it from the main body. (f)
4. Empty the receptacle of callused skin and use the cleaning brush to finish. (g)
5. To re-assemble the device, first connect the receptacle to the body of the device, verifying the direction of the magnet. (h)
6. Place the abrasive disk on the receptacle, making sure to properly secure it. (i)
7. Align the circular bulge of the main body with that of the safety ring and then turn the ring counter-clockwise. (j & k)
8. Put the cover back on if you are not using the device. (l)

MAINTENANCE AND STORAGE:

MAINTENANCE

- Clean all of the parts after each use as described above and use a moist cloth for the exterior of the device.
- Do not use detergent or solvent.

STORAGE

After cleaning and drying the device, store it in a dry, cool place out of the reach of children.

ADVICE REGARDING THE DISPOSAL OF WASTE

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control. LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

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