- 1: Temperatursensor
- 2: Display
- 3: On (Ein)/Start-/Stop-/Speicher-/Off (Aus)-Taste

Temperaturanzeige

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
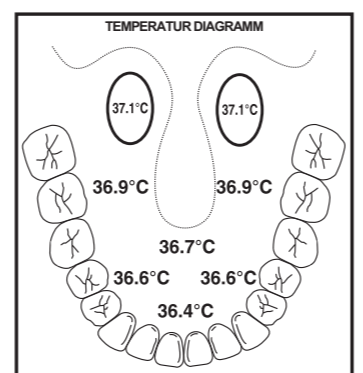
2. Place the thermometer in the chosen body part to measure your temperature.
4. When you hear beeps, your temperature has been recorded.
 - if your body temperature is normal, you will hear 10 beeps
 - if your body temperature is above normal (≥37.8°C), you will hear 10x3 beeps.
5. Press the ON/OFF key to switch off the thermometer. (The thermometer will switch off automatically after approx. 10 minutes).

*: If the room temperature is higher than 32.0°C, the room temperature will be shown instead of 'Lo'.

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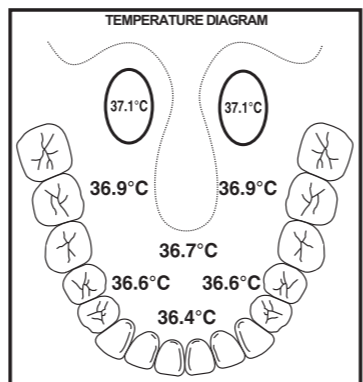Measurement under the arm
 Make sure the armpit is dry and place the thermometer under the arm. The sensor must touch the skin. Keep your arm crossed in front of your chest. A temperature taken under the tongue is approx. 1.0°C lower than when taken rectally.

Rectal measurement
 Insert the thermometer sensor gently (maximum 1.25 cm).

REPLACING THE BATTERY

The battery must be changed when the symbol 'LO' appears in the display. Battery type: 1.55V, LR41 alkaline or SR41 silver oxide type or similar.

1. Carefully remove the battery cover from the thermometer.
2. Use a pointed (non-metal) object to remove the battery. Insert the new battery with the + side upward.
3. Slide the battery cover back in place.

Do not throw away empty or leftover batteries, but instead hand them in to your local depot for Household Chemical Waste.

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ADVICE AND WARNINGS

- Keep the thermometer and battery out of the reach of children.
- Do not talk or walk around while taking your temperature.
- Store the thermometer in the special case.
- Do not expose the thermometer to direct sunlight, dust or moisture.
- Do not dismantle the thermometer.

STORAGE AND CLEANING

Clean the thermometer before and after use, using a soft cloth and alcohol, or cold soapy water. The thermometer may be briefly immersed in water or alcohol (do not soak or sterilise). Store the thermometer in the special case.

CE

The declaration of conformity is available at:
<http://DOC.hesdo.com/BC-13-DOC.pdf>

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
WARRANTY

For the Alecto BC-13 thermometer you have a warranty of 24 months from the date of purchase. During this period, we guarantee the free repair of defects caused by material and workmanship errors. All this subjected to the final assessment of the importer.

HOW TO HANDLE: If you notice any defect, first refer to the user's manual. Please contact the supplier of this thermometer or contact customer service of Alecto via the Internet www.alecto.nl

THE WARRANTY BECOMES NULL AND VOID: In case of improper use, incorrect connections, leaking and/or incorrectly installed batteries, use of unoriginal parts or accessories, negligence and in case of defects caused by moisture, fire, flooding, lightning and natural disasters.

In case of unauthorised modifications and/or repairs performed by third parties. In case of incorrect transportation of the device without adequate packaging and when the device is not accompanied by this warranty card and proof of purchase. All further liability, especially concerning any consequential damages, is excluded.

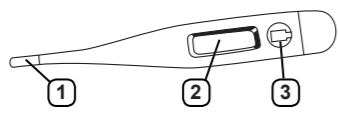- 1: temperature sensor
- 2: display
- 3: on/start/stop/memory/off button

temperature display

temperature display memory

battery empty

INTRODUCTION

This Alecto thermometer is a digital thermometer developed to take your temperature as simply and accurately as possible. In order to get the most reliable results, these instructions must be followed accurately. For the most accurate temperature reading, we recommend using the thermometer rectally. Although the generally accepted 'normal' temperature is 37.0°C, body temperature can vary from 36.1°C to 37.2°C.

Fluctuations in temperature can be due to activities such as sports, smoking, eating and drinking, and also depend

on the part of the body where the temperature is taken. When measured orally, the body temperature is generally 0.5°C lower than when measured rectally, while a temperature measured under the armpit is 1.0°C lower. Your temperature is lower in the morning than in the evening.

SPECIFICATIONS

Temperature range: 32.0°C - 42.9°C
 Low temperature display: temperature ≤ 32°C, display: Lo°C
 High temperature display: temperature ≥ 42.9°C, display: HI°C
 Display resolution: 0.1°C
 Accuracy: 35.5°C - 42°C: ± 0.1°C
 ≤ 35.5°C ≥ 42°C: ± 0.2°C
 Normal temperature ≤ 37.8°C: 10 beeps
 Fever alarm ≥ 37.8°C: 10x3 beeps
 Auto-OFF: 10 minutes ± 2 minutes
 Battery: LR41, 192 or SR41
 Power consumption: 0,15 mW
 Battery life: 2 years
 Dimensions: 126 x 19 x 12,5mm
 Weight: ±10 gr (batt incl.)
 Conditions for use:
 Temperature: 10°C - 35°C
 Humidity: 25% - 80% RH
 Storage conditions:
 Temperature: -25°C - 55°C
 Humidity: 25% - 80% RH

MEASURING YOUR TEMPERATURE

1. Press the ON/OFF key once to activate the thermometer. You will hear a beep.
2. The display shows in succession: '188,8°C' the last measured temperature '37,0' (for internal calibration) 'LO' and flashing '°C' (*)

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Hangzhou Universal Electronic Co., Ltd
Prolinx GmbH